

---

---

# ALCOSCAN AL3500

*Coin-Operated Breath Alcohol Tester*



## OPERATING MANUAL



# CONTENTS

1. Introduction
2. Parts and Functions/Display Panel
3. Device Operation
4. Warnings/Precautions
5. Maintenance
6. Troubleshooting



# INTRODUCTION

The AlcoScan AL3500 is a coin-operated breathalyzer (breath alcohol tester) suitable for use in bars, restaurants and other commercial locations.

Interactive voice prompts (integrated ARS/ automatic voice response system) and bright LED indicators on the front of the device provide a quick and easy way for users to check their alcohol levels (in %BAC).

The AL3500 features advanced oxide semiconductor sensing technology, as well as patent-pending Pre-Calibrated Replaceable Sensor Modules that can be used in place of calibration.


**IMPORTANT: The only legally valid measure of %BAC (percent Blood Alcohol Concentration) is obtained by blood testing. This device is for screening purposes only and should be used as a guide for self-monitoring.**

# PARTS AND FUNCTIONS



- 01 Mouthpiece Slot ("Bullseye")
- 02 Coin Slot
- 03 Display Panel
- 04 Straw Tray
- 05 Fan
- 06 Speaker
- 07 Power Socket

## DISPLAY PANEL

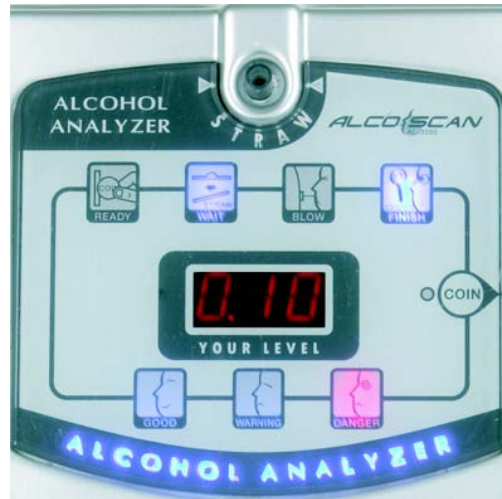

01 "Ready" Indicator 02 "Wait" Indicator 03  
"Blow" Indicator 04 "Finish" Indicator 05 LED  
Alcohol Level Display (%BAC) 06 "Good"  
Indicator (%BAC < 0.02) 07 "Warning"  
Indicator (%BAC from 0.02 ~ 0.05) 08 "Danger"  
Indicator (%BAC > 0.05)

# DEVICE OPERATION


**Step 1:** Insert the appropriate number of coins into the coin slot.

Voice prompt: **“Thank you. Please take a straw from the tray and wait until the machine counts down to ‘0’ in the display.”**

**Step 2:** Take a straw from the tray and insert it into the mouthpiece slot (“bullseye”). When the device indicates “Blow”, blow into the mouthpiece steadily until the device indicates that sampling is completed. If a successful sample is not taken, the device will prompt you to wait, then blow another breath sample until sampling is successful.

Voice prompt: **“Insert the straw securely into the ‘bullseye’ until the tone ends. Three, two, one <CHIME SOUND>. Sampling is completed.”**

Voice prompt (if successful sample is not taken, usually due to weak breath): **“The device does not sense a proper breath sample. Please wait and try again.”**



**Step 3:** After a successful sample is taken, the device will analyze the sample and display the alcohol level as %BAC (% Blood Alcohol Concentration). In addition to displaying the %BAC, an LED indicator will also light up, as described below.

**LED Indicator "GOOD" (%BAC less than 0.02)**

Voice prompt: **"The device detects very little or no alcohol. However, if you have had ANY amount of alcohol, it may not be safe to drive."**

**LED Indicator "WARNING" (%BAC between 0.02 and 0.05)**

Voice prompt: **"The device detects a moderate level of alcohol. Please wait at least an hour or more and try another test. It may not be safe to drive."**

**LED Indicator "DANGER" (%BAC greater than 0.05)**

Voice prompt: **"The device detects a high level of alcohol. It is NOT safe to drive."**



**IMPORTANT NOTE:**

**If the user does not blow into the straw after coins are inserted and the countdown runs to "0", a ten-second timer will appear in the display. If a breath sample is not given and the timer runs out, the test can no longer be taken. Coins will not be refunded.**



## WARNINGS/PRECAUTIONS

**Make sure you carefully read this section before using the AL3500.**

Please wait at least 20 to 30 minutes after drinking, eating or smoking before performing a breath test with the AL3500. Liquid alcohol, cigarette smoke, excess saliva and other contaminants in the mouth can greatly affect the test results.

Do NOT use the AL3500 as a tool to drink and drive.

The only legally valid measure of %BAC (percent Blood Alcohol Concentration) is obtained by blood testing. The AL3500 is for screening purposes only and should be used as a guide for self-monitoring.

# MAINTENANCE

The AL3500 offers powerful features and performance as well as convenient, user-friendly configuration. There are many customizable functions available.

In order to configure the AL3500, open the rear door of the device and locate the three buttons along the left side of the circuit board.

The buttons are labeled **SET1**, **SET2** and **TEST** from top to bottom.

**NOTE:** The **TEST** button automatically starts a breath test on the AL3500, equivalent to inserting coins.



**SET1 (SW1)**



**SET2 (SW2)**



**TEST (SW3)**



## CUMULATIVE TEST COUNTER

In order to check the total (cumulative) number of tests performed on the AL3500, press **SET1** for one second. The LED display on the front of the device will show the number of tests.

Press and hold **SET1** for two seconds to reset this test counter to zero.

If no commands are entered for several seconds, the device and display will return to normal operating mode.

## COIN ACCEPTOR CONFIGURATION

In order to adjust the number of coins required to take a breath test on the AL3500, press **SET2** for one second. The LED display on the front of the device will show a number from "**C01**" to "**C10**", depending on the current setting for the number of coins (one through ten).

Tap the **SET2** button to change the number of coins.

If no commands are entered for several seconds, the device and display will return to normal operating mode.

## **ADVANCED CONFIGURATION MODE (Buzzer, Sound, Lo/Hi Levels)**

If you want to adjust the buzzer, sound and Lo/Hi Levels, you must access the Advanced Configuration Mode. This is done by simultaneously pressing and holding both **SET1** and **SET2** for two seconds.

While you are in Advanced Configuration Mode, tapping **SET1** will cycle through the options of Buzzer, Sound, Lo and Hi. **SET2** will change the setting for each option.

### **Buzzer**

Tap **SET1** until you see “**buz**” in the LED display. Tap **SET2** to turn the buzzer (audio beeps) on or off.

### **Sound**

Tap **SET1** until you see “**Snd**” in the LED display. Tap **SET2** to turn the interactive voice prompts on or off.



### **Alcohol Level - Low Limit**

This setting adjusts the detected alcohol level below which the AL3500 will light the “Good” Indicator after a test is performed.

The default setting is **0.02**. **We do not recommend changing this setting.**

Tap **SET1** until you see “**Lo**” in the LED display. Tap **SET2** to adjust the limit from **0.00** to **0.04**.

### **Alcohol Level - High Limit**

This setting adjusts the detected alcohol level above which the AL3500 will light the “Danger” Indicator after a test is performed.

The default setting is **0.05**. **We do not recommend changing this setting.**

Tap **SET1** until you see “**Hi**” in the LED display. Tap **SET2** to adjust the limit from **0.05** to **0.10**.

## REPLACING SENSORS

The AL3500 is equipped with Pre-Calibrated Replaceable Sensor Modules. Although every consumer breathalyzer requires periodic re-calibration, the AL3500 allows you to simply discard an old or worn sensor and replace it with a brand new factory calibrated sensor module.

Sensor Replacement is recommended every 250~300 tests, or at least once a year (whichever is sooner).

Sensor modules can also be replaced if there is defect or malfunction in the sensor.

### Sensor Module Replacement Instructions

**Step 1:** Power Off the AL3500 (**VERY IMPORTANT - The new sensor may be severely damaged if power is on**).

**Step 2:** Remove the old sensor module from the AL3500 circuit board and replace it with a new sensor module.

**Step 3:** Make sure the sensor module is seated securely, and power on the AL3500.



# TROUBLESHOOTING

## ERROR CODES

There are two error codes that may be shown in the LED display in cases of malfunction or other problems.

**Er3** - This code may appear if an improper coin is inserted into the coin slot, causing it to get stuck inside the coin acceptor. Try turning the coin release lever or examining the coin acceptor by opening the rear door of the device.

**SEn** - This code may appear if the sensor module is not installed properly, or if there is a problem or malfunction in the sensor. You will most likely need to replace the sensor module.

## COMMON PROBLEMS

### 1. The AL3500 does not power on

Make sure that the power cable is securely attached to both the wall socket as well as the power socket at the bottom of the AL3500. Also, check inside the AL3500 to make sure that the wiring from the power socket is properly connected (red/black wires at the bottom).

## **2. Coins are inserted, but the test doesn't activate.**

Make sure you are using the proper type and number of coins. Also, check your coin acceptor setting, as described on Page 10 in this manual.

## **3. When you blow into the machine during a test, there is no response from the machine**

See if the LED display shows error code "**SEn**". If so, check the error codes on Page 14 in this manual. Another common problem is that the sensor module is not seated correctly. Open the rear door of the AL3500 and make sure that the sensor module is securely attached to the black breath tube at the top of the circuit board.

## **4. There are no interactive voice prompts and/or audio**

The voice prompt and audio options can be adjusted in Advanced Configuration Mode. See Page 11 in this manual.



**5. After a breath test, the AL3500 does not return to normal operating mode.**

After each test, the AL3500 displays the results for approximately ten seconds. In order to return to normal operating mode (for a new test) the alcohol residue from the current test must be removed completely. This is an automatic process, but may take longer if there is a high alcohol reading (up to one minute). If this takes much longer than one minute, there may be a problem with the sensor assembly. Contact technical support using the information below to resolve your issue.

**6. I have a problem that is not covered in this Troubleshooting section.**

Please contact technical support:

**866-216-8700 Toll Free**  
**818-701-9200 California**  
**818-701-9220 FAX**

Distributed By:

**Advanced Safety Devices  
9140 Jordan Avenue  
Chatsworth, CA 91311**

**Tel: 866-216-8700 US  
Tel: 818-701-9200 Intl.  
Fax: 818-701-9220**

**[sales@safety-devices.com](mailto:sales@safety-devices.com)**

**[www.safety-devices.com](http://www.safety-devices.com)**